

Enhanced convenience:
UAW ClearShield



Less splashes and mist during wound debridement:

UAW ClearShield

Soft, flexible silicone cap enhancing convenience during UAW through:

Special design:

Splashes and mist adhere to the inside of the cap, enhancing convenience during UAW and reducing the amount of needed wipe-disinfection. Shape enables sufficient view on wound from any angle.

Unique coating:

Upon activation with fluid, the special coating becomes transparent and ensures clear view on wound without fogging.

Soft silicone material:

Can easily be moved during application. High-quality silicone has low risk of allergies and is very skin-friendly.



UAW FlexLock:

Safely connects UAW ClearShield with UAW Instrument. Flexible, accordion-like composition allows for utmost application comfort and tactile feedback.




Ultrasonic-Assisted Wound Debridement (UAW) with SONOCA 185 from Söring enables safe and selective wound debridement and biofilm disruption, whereas healthy tissue is protected, making it very beneficial for debridement and cleansing of chronic wounds and also for preparation of skin grafts.

UAW ClearShield enhances convenience during Wound Debridement as it effectively reduces splashes and mist during treatment. It is compatible with Söring UAW instruments through an adapter. Soft and safe UAW ClearShield can easily be moved during treatment and can be used with one or two hands in order to suit a variety of wounds. The special shape makes UAW ClearShield ideal for venous leg ulcers and pressure ulcers.




Less splashes and mist with UAW ClearShield

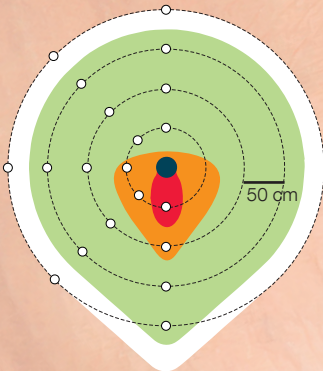

Spread of splashes and mist without UAW ClearShield ¹

- Sonotrode tip
- Measuring point
- 0 - 1.25 particle/mm²
- 1.25 - 7.48 particle/mm²
- 7.48 - max. particle/mm²

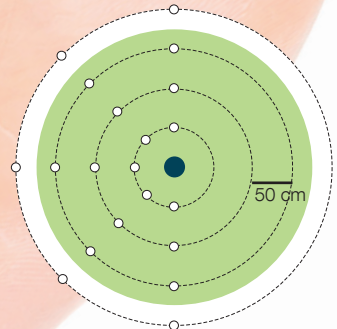

Spread of splashes and mist with UAW ClearShield ¹

¹ Söring GmbH laboratory data

Product overview

Article no.

UAW ClearShield, 10 pcs per folding box, single-use

200N0104-10

UAW ClearShield connection set for UAW instruments

97-102, 97-103, 97-104 reusable*, contains:

615S0056

- UAW FlexLock

615S0058**

- UAW FlexLock Adapter

615S0059**




* after reprocessing

** not separately available

UAW ClearShield
(single-use)

UAW FlexLock
(reusable)

UAW FlexLock Adapter
(reusable)

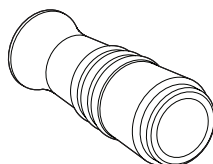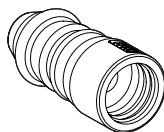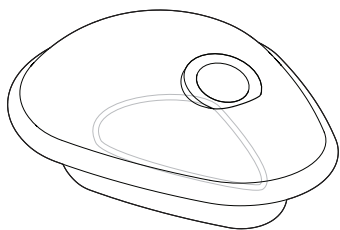

About Söring

Söring has been developing innovative medical technology “Made in Germany” for various surgical disciplines for over 30 years. The family-run business, founded and still located in Quickborn (near Hamburg), utilizes extensive experience in ultrasonic technology, a dedicated in-house R&D Team and close clinical collaborations for new solutions to support clinicians in their daily clinical practice.

→ Söring GmbH | Justus-von-Liebig-Ring 2 | 25451 Quickborn | Germany
Phone: +49 4106-6100-0 | E-mail: info@soering.com

Further information can be found at: www.soering.com

Specifications, design and accessories of the shown products might have changed in the meantime.

The products shown might not be available in your country. Please contact the local Söring partner for further information.