

ENTERPRISE 5000 HOSPITAL BED



Inspired design can be refreshingly simple!

The new Enterprise® 5000 bed is the latest general use hospital bed in the Enterprise range. The new design may offer both caregiver and patient benefits in the key areas of safety, clinical efficacy and ease of use.

A well-designed hospital bed can offer significant advantages. These include enhanced safety, increased patient

independence and improved efficiency. A hospital bed needs to be as user friendly as possible, enabling the caregiver time to concentrate on their patient.

Enterprise 5000 bed is easy and intuitive to use, thanks to clear, simple icons and controls designed to make the operation of bed functions as straightforward as possible.




...with people in mind

www.ArjoHuntleigh.com


KEY FEATURES

2 3 Step Bed Extension

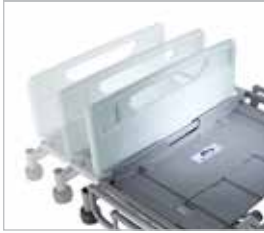

Different bed lengths are now easily configurable in a single sub-frame, offering a simple to use, multi-length hospital bed for greater flexibility and choice.

3 Underbed Light

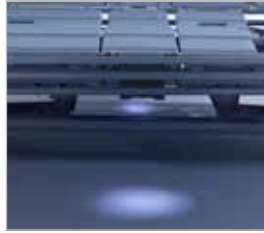

Optional underbed lighting fitted under both sides of the bed helps to improve patient safety and security, by illuminating the bed exit area.

4 Ease of Access



The open architecture design allows ease of access to help improve decontamination procedures. Deck sheets are lightweight, smooth and easily removed.

1 Curved Mattress Deck



The patented *Bio-Contour* profiling system blends a unique combination of bed movements with a curved mattress deck.



Enterprise 5000 bed also has an extensive offering of accessory items, which both complements the bed and provides additional benefits for patients and carers. Please contact your local ArjoHuntleigh representative for further details.

5 Pit Stop Maintenance



Designed for 'Pit Stop Maintenance'- open accessibility means changes such as actuator, control box and battery can be carried out in situ within a few minutes by one person.

6 Attendant Control Panel



All new One-Button-Push control panel with updated graphics and layout for intuitive ease of use. Features auto-CPR, lock-outs and cardiac chair position.

7 Patient Handset



A choice of hand controls is available to provide independent comfort adjustment for the patient.



8 30 Degree Backrest Angle



The backrest will pause when raised to an angle of 30 degrees, ensuring easy and positive back-rest adjustment time after time.

9 Low Height



An exceptional low height of 12.5ins (320mm) is achieved via One-Button-Push.

PRODUCT SPECIFICATIONS

FEATURES	
Electrically operated backrest, knee-break and height adjustment	•
Electrically operated Trendelenberg / Reverse Trendelenberg tilt	±12°
Auto-CPR and Auto-contour button	•
Battery back-up	•
Bio-Contour profiling system	•
One button chair position	•
Dual sided CPR levers	•
30° backrest angle pause	•
Angle indicators for backrest	•
Built-in bed extension with mattress platform extension	•
Removable head and foot panels with choice of colour infills	•
Drainage bag rails	•
Removable plastic deck sheets	•
Linked braking / steering system with four braking castors	•
125mm single wheel castors	(•)
150mm single wheel castors	(•)
150mm twin wheel castors	(•)
Additional brake pedals at head end	(•)
Full width brake bar (foot end only)	(•)
DIN rails	(•)
Under-bed light	(•)
Folding safety sides	(•)
Pull-out linen shelf	(•)

• Standard (•) Optional

SPECIFICATIONS	
Overall standard length	2300mm
Overall width	1030mm
Mattress length standard	2020mm
Mattress width	880mm
Overall extended length	2420mm
Mattress extended length	2140mm
Overall length short	2190mm
Mattress length short	1910mm
Top height (with 125mm castors)	760mm
Low height (with 125mm castors)	320mm
Backrest length	850mm
Seat length	215mm
Calf length	560mm
Thigh length	365mm
Safe working load/Maximum patient weight	250kg/185kg
Head down tilt	12°
Foot down tilt	12°
Backrest angle	62°
Backrest pause	30°
Thigh angle	20°
Calf angle	16°
Ingress protection	IPX4
Shock protection	Class 1, Type B
Power in	1.25A at 230V AC 50 / 60Hz
Power in	2A at 120V AC 50 / 60Hz
Safety standard	IEC60601-2-52

ArjoHuntleigh is a branch of Arjo Ltd Med. AB. Only ArjoHuntleigh designed parts, which are designed specifically for the purpose, should be used on the equipment and products supplied by ArjoHuntleigh. As our policy is one of continuous development we reserve the right to modify designs and specifications without prior notice.

ARJOHUNTLEIGH
GETINGE GROUP

GETINGE GROUP is a leading global provider of products and systems that contribute to quality enhancement and cost efficiency within healthcare and life sciences. We operate under the three brands of ArjoHuntleigh, GETINGE and MAQUET. ArjoHuntleigh focuses on patient mobility and wound management solutions. GETINGE provides solutions for infection control within healthcare and contamination prevention within life sciences. MAQUET specializes in solutions, therapies and products for surgical interventions and intensive care.

Only ArjoHuntleigh designed parts, which are designed specifically for the purpose, should be used on the equipment and products supplied by ArjoHuntleigh. As our policy is one of continuous development we reserve the right to modify designs and specifications without prior notice.

® and ™ are trademarks belonging to the ArjoHuntleigh group of companies.
© ArjoHuntleigh, 2013.

www.ArjoHuntleigh.com



ArjoHuntleigh AB
Hans Michelsensgatan 10
SE-211 20 MALMÖ
SWEDEN
Phone: +46 10 335 45 00